



NEW!

AeroFlex™ For

Printed Resins • PMMA • Pekkton • Peek • Clear Aligners • Acrylics • Thermo Materials

AeroFlex™ has a reinforced binder that provides positional stability allowing for controlled area contouring, deburring and flash removal. AeroFlex™ quickly eliminates scratches and leaves smooth clearly defined surfaces. In short order your restoration will be ready for a final high shine polish.

AeroFlex™



Mastering Green State Zirconia

The GS Green State Contour™ System provides a variety of instruments that offer a safe and efficient way to remove sprue connectors and contour your pre sintered zirconia restorations.



GS Kit



GS Contour Diamond Blade

Our diamond coated blade was designed to open and contour interproximal areas. The blade can gently access tight areas without causing any damage. Blades are available in fine and ultra fine.



GSBK Diamond Blade Kit

- 1 Handle
- 5 Blades Ultra Fine
- 5 Blades Fine

ZHP Diamonds

Our 6 ZHP hand piece diamonds have a unique diamond shape and particle size. They provide a safe, gentle way to make subtle anatomical adjustments in the Green State.





Knife Edge
26x2mm



Golden Eagle
For Feldspathic Ceramics

Our Golden Eagles are infused with a high percentage of diamond and contour feldspathic ceramic beautifully. We also removed the yellow centers to give you more usable disc.

Golden Eagle Kit

26mm: Extra coarse, coarse, medium, fine-high shine
17mm: coarse, medium, fine
Points: coarse, fine



ITEM GEK1

Flat Edged Wheels
17x3mm



Large Point
5.5x18mm




Small Point
3.3x7.5mm



Diamond Discs For Zirconia

While a traditional diamond disc works well with high fusing feldspathic ceramic, you have probably noticed that it's less effective and does not last as long when used with zirconia or lithium disilicate. Our ZR discs have a special diamond shape, grit, and tighter particle placement. They are more effective, last longer and stay cooler. Our 38mm and 30mm discs are recommended for cutting pressed sprues. Available Disc Diameters: 08,10,14,16,19,22,30,38,40,45

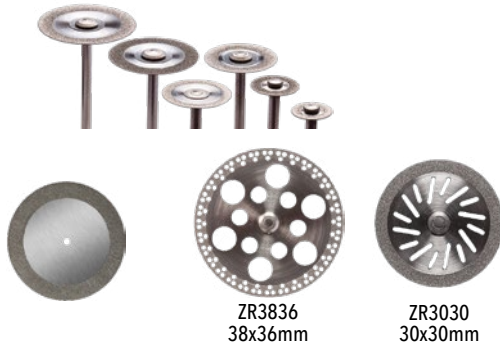

ZRHP Diamonds
For All Ceramic Restorations

Our ZR Hand Piece diamonds provide a safe alternative to using FG diamonds when irrigation is not available. They are the perfect solution for final, minor contouring of all ceramic restorations.




Chairside
Contour...Adjust...Polish...Or Remove

Zirconia, Lithium Disilicate, Nano Hybrid Composition



RAFG Kit

- 2 RA RedBerry Coarse
- 2 RA BlueBerry Medium
- 4 RA GoldenBerry Fine High Shine
- 7 FG For Adjusting
- 2 FG For Removal

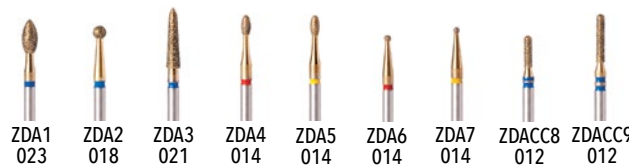



Finale™ Final Contouring System

Finale's uniquely shaped fine diamond particles allow for gentle final contouring without heat transfer, and are particularly well suited for contouring thin margins. Finale™ points permit access to difficult to reach areas. Finale™ instruments are available in fine and medium grits; use Finale after using your RedBerry® or BlueBerry® instruments.




F Kit



Friction Grip Diamonds For Zirconia & Lithium Disilicate Restorations

Our ZR FG diamonds have been layered with uniquely shaped particles that make them highly effective and safe to use with your all-ceramic restorations. Only minimal hand pressure is required to obtain maximum results. Our FG diamonds are long lasting, do not leave black marks, and are available in a variety of shapes and grits. Our ZR FG diamonds should always be used with irrigation to prevent heat transfer that can create micro fractures in your restorations. We can also produce custom shapes in a variety of grits to meet your specific needs.

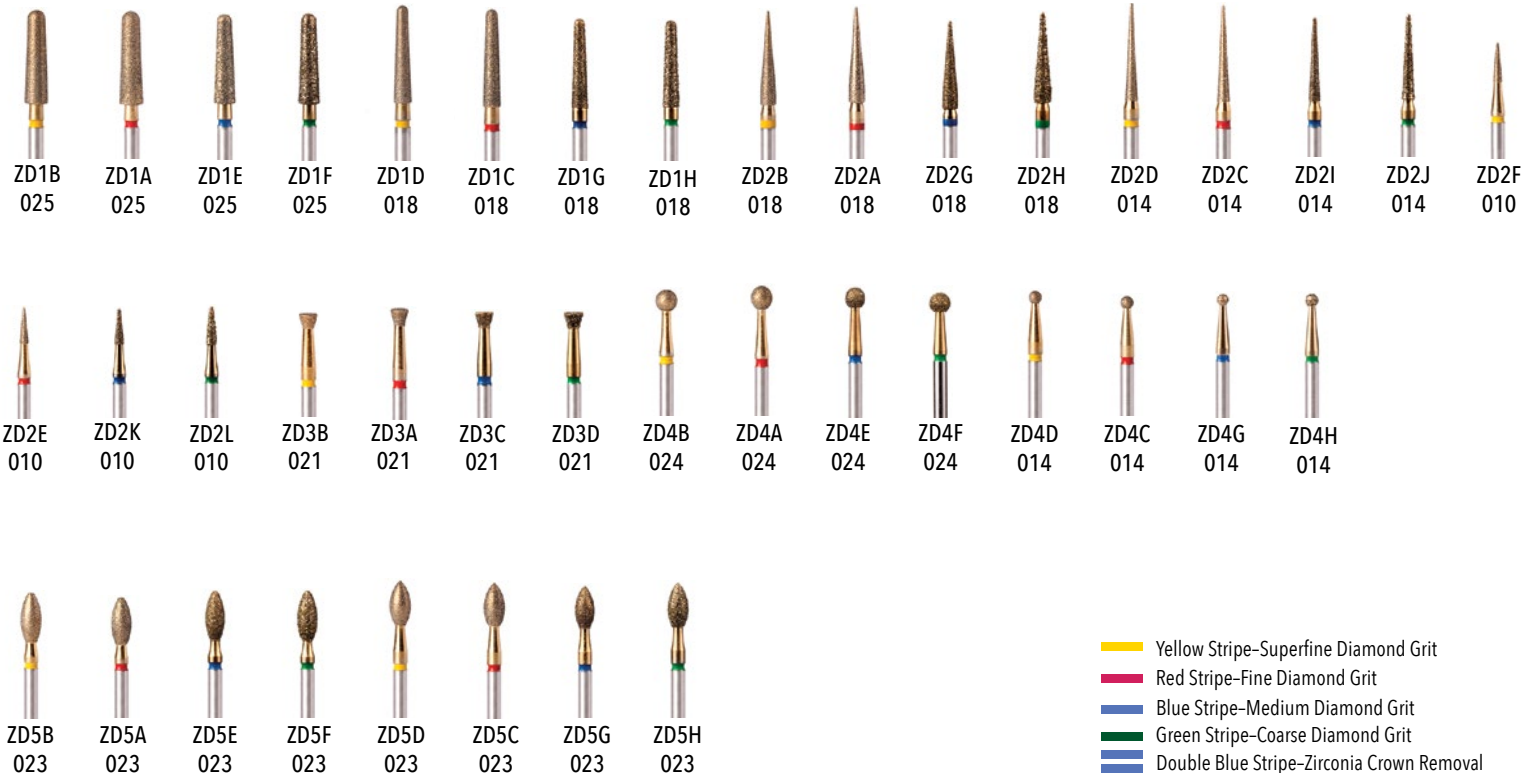

NEW!



When a greater amount of material needs to be removed, our new Diacool ICE instruments will quickly and easily get the job done. Unique diamond particle shape and size, combined with a new binder, ensure zirconia, lithium disilicate or feldspathic restorations remain cool, thereby reducing the possibility of creating micro fractures while grinding. With only minimal pressure applied, restorations are quickly and safely reduced and

shaped. Some ceramists prefer using the coarse version of Diacool ICE; just slightly more abrasive, yet they also remain cool and effective. The large B406 flat edged wheel is an excellent choice for pressed lithium disilicate sprue connector reduction. For final contouring we recommend using our "Berrys", or "Finale" for adjusting thin margins.

Smoother, Gentle Material Removal • Minimal Heat Transfer • Increased Instrument Life



Medium



Coarse



DCI Kit

The "Berrys"

Our Berry instruments were specifically developed for use with zirconia and lithium disilicate restorative materials. They are infused with a high concentration of uniquely shaped diamond particles that allows for more effective contouring and shaping. The "Berrys" reduce heat transfer, prevent microscopic surface scratches, and leave a smooth pre-polish surface. This smooth surface eliminates the capillary effect that causes stains to be drawn to unwanted areas; that's industry unique, and we are certain you have never experienced diamond infused instruments that perform like our "Berrys".




UltraBerry® Coarse Plus, Diamond Infused Preferred for zirconia contouring




BK Kit



GoldenBerry® Super Fine, Diamond Infused. Places a fast super high shine on zirconia or lithium disilicate restorations

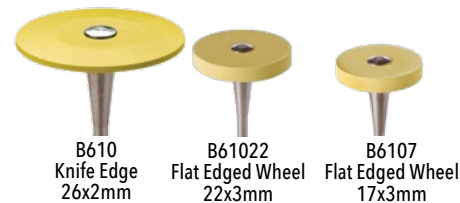

RedBerry® Coarse, Diamond Infused; preferred for zirconia contouring

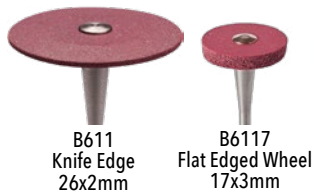

NEW!



BlueBerry® Medium, Diamond Infused; preferred for lithium disilicate contouring, or as a second step when contouring zirconia



NEW!



NEW!



NEW!

Occlusal Contouring and Polishing

Wagner Diamond Infused Berry Spirals provide a flexing action that allows fingers to get into difficult interproximal and occlusal areas.

- RedBerry® & BlueBerry® Spirals are used to make final adjustments and refine anatomy
- GoldenBerry® Spiral will create a high shine polish

